

## GMP PE/Cyanine5 anti-human CD33 Antibody

<b>Catalog# / Size</b>	260366 / 100 tests
<b>Clone</b>	WM53
<b>Workshop</b>	IV M-505
<b>Other Names</b>	Siglec-3, gp67, p67
<b>Isotype</b>	Mouse IgG1, κ
<b>Description</b>	CD33 is a 67 kD type I transmembrane glycoprotein also known as Siglec-3, gp67, and p67. It is a sialoadhesion immunoglobulin superfamily member expressed on myeloid progenitors, monocytes, granulocytes, dendritic cells and mast cells. CD33 is absent on normal platelets, lymphocytes, erythrocytes and hematopoietic stem cells. CD33 functions as a sialic acid-dependent cell adhesion molecule with carbohydrate/lectin binding activity.

### Product Details

<b>Reactivity</b>	Human
<b>Reported Reactivity</b>	Chimpanzee
<b>Antibody Type</b>	Monoclonal
<b>Host Species</b>	Mouse
<b>Immunogen</b>	Human myeloid leukaemia cells.
<b>Formulation</b>	Phosphate-buffered solution, pH 7.2, containing True-Stain Monocyte Blocker™, 0.09% sodium azide and 0.2% (w/v) BSA (origin USA) and a stabilizer.
<b>Preparation</b>	The antibody was purified by affinity chromatography and conjugated with PE/Cyanine5 under optimal conditions.
<b>Concentration</b>	100 µg/mL
<b>Storage &amp; Handling</b>	The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. <b>Do not freeze.</b>
<b>Application</b>	<a href="#">FC - Quality tested</a>
<b>Recommended Usage</b>	Each lot of this antibody is quality control tested by <a href="#">immunofluorescent staining with flow cytometric analysis</a> . For flow cytometric staining, the suggested use of this reagent is 5 µL per million cells in 100 µL staining volume or 5 µL per 100 µL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.
<b>Excitation Laser</b>	Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm)
<b>Application Notes</b>	Additional reported applications (for the relevant formats) include: immunoprecipitation, Western blotting <sup>3</sup> , induction of cytokine production <sup>3</sup> , and immunofluorescence <sup>4</sup> . The Ultra-LEAF™ purified antibody (Endotoxin < 0.01 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for functional assays (Cat. Nos. 303437 & 303438).
<b>Application References</b>	<ol style="list-style-type: none"> <li>Knapp W, <i>et al.</i> 1989. Leucocyte Typing IV. Oxford University Press. New York.</li> <li>Favaloro E, <i>et al.</i> 1988. <i>Br. J. Haematol.</i> 69:163.</li> <li>Garnache-Ottou F, <i>et al.</i> 2005. <i>Blood</i> 105:1256. (WB)</li> <li>Pérez-Oliva AB, <i>et al.</i> 2011. <i>Glycobiology.</i> 21:757. (epitope, FC, IF)</li> </ol>
<b>(PubMed link indicates BioLegend citation)</b>	
<b>Disclaimer</b>	<b>GMP RUO Flow Cytometry Antibodies.</b> BioLegend GMP RUO fluorophore conjugated antibodies are manufactured in a dedicated GMP facility and compliant with ISO 13485:2016. For research use only. Not for use in diagnostic or therapeutic procedures. Our processes include: <ul style="list-style-type: none"> <li>• Batch-to-batch consistency</li> <li>• Material traceability</li> </ul>

- Documented procedures
- Documented employee training
- Equipment maintenance and monitoring records
- Lot-specific certificates of analysis
- Quality audits per ISO 13485:2016
- QA review of released products

## Antigen Details

<b>Structure</b>	Ig superfamily, sialoadhesins, type I glycoprotein, 67 kD
<b>Distribution</b>	Myeloid progenitor, monocytes, granulocytes, dendritic cells, mast cells
<b>Function</b>	Adhesion and lectin activity
<b>Ligand/Receptor</b>	Sugar chains containing sialic acid
<b>Cell Type</b>	B cells, Dendritic cells, Granulocytes, Hematopoietic stem and progenitors, Mast cells, Monocytes, Neutrophils
<b>Biology Area</b>	Cell Biology, Immunology, Neuroinflammation, Neuroscience
<b>Molecular Family</b>	CD Molecules, Siglec Molecules
<b>Antigen References</b>	1. Favaloro E, <i>et al.</i> 1988. <i>Br. J. Haematol.</i> 69:163. 2. Freeman S, <i>et al.</i> 1995. <i>Blood</i> 85:2005.
<b>Gene ID</b>	<a href="#">945</a>


## Related Protocols

[Cell Surface Flow Cytometry Staining Protocol](#)

## Other Formats

APC anti-human CD33, PE anti-human CD33, PE/Cyanine5 anti-human CD33, Purified anti-human CD33, PerCP/Cyanine5.5 anti-human CD33, Brilliant Violet 711™ anti-human CD33, Brilliant Violet 421™ anti-human CD33, Brilliant Violet 570™ anti-human CD33, Purified anti-human CD33 (Maxpar® Ready), Brilliant Violet 510™ anti-human CD33, Biotin anti-human CD33, Brilliant Violet 785™ anti-human CD33, Brilliant Violet 650™ anti-human CD33, PE/Cyanine7 anti-human CD33, PE/Dazzle™ 594 anti-human CD33, Alexa Fluor® 700 anti-human CD33, Ultra-LEAF™ Purified anti-human CD33, APC/Fire™ 750 anti-human CD33, APC/Cyanine7 anti-human CD33, GMP PerCP/Cyanine5.5 anti-human CD33, Alexa Fluor® 647 anti-human CD33, GMP APC anti-human CD33, GMP PE anti-human CD33, PE/Fire™ 810 anti-human CD33, Spark Blue™ 515 anti-human CD33, PE/Fire™ 700 anti-human CD33

## Product Data

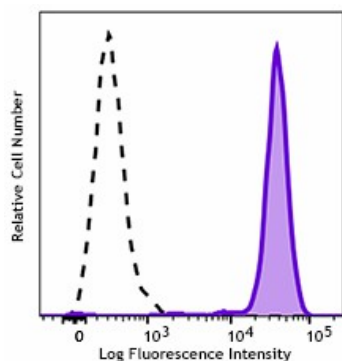

Typical results from human peripheral blood monocytes stained either with WM53 PE/Cyanine5 used at 5  $\mu$ L/test (filled histogram) or with an isotype control (open histogram).

\*These products may be covered by one or more Limited Use Label Licenses (see the BioLegend Catalog or our website, [www.biolegend.com/ordering#license](http://www.biolegend.com/ordering#license)). BioLegend products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products, reverse engineer functionally similar materials, or to provide a service to third parties without written approval of BioLegend. By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses.

8999 BioLegend Way, San Diego, CA 92121 [www.biolegend.com](http://www.biolegend.com)  
Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587